

# HOME CARE LINE

sapimed

CE

DEPLPRODING REV02 09 2018



**Dilatan Plus** is a cryo-thermal anal dilator available in 5 different sizes only to be used when advised by a physician. Three boxes are available (18/20 mm, 23/27 mm, 30 mm), offering patients the most suitable size for their particular problem and individual anatomy. Dilatan Plus consists of an inert plastic shell with a large disk at the base to allow the patient to grip the device easily. Each sealed dilator contains a thermal jelly which, when heated or cooled, maintains the dilator temperature during use.

REFs : • A.1001 • A.1002 • A.1030



The **DILASTOM®** stomal dilator is a medical device specifically designed so that patients may carry out home treatment for cutaneous and fascial colostomy stenosis with complete safety and convenience. Surgical operations to correct colostomy stenosis may be avoided by protracted use of DILASTOM® for long periods, without producing any undesirable side effects.

REFs : • A.2002

**DILAGENT** is a **soft silicone** anal dilator. It is indicated for the treatment of anorectal diseases caused by a hypertonic sphincter, namely anal fissures, haemorrhoids and painful spasms after surgical treatment of the anorectal segment. It is also effectively used in cases of postsurgical stenosis of the anal canal. It is successfully used to treat Vaginism and Dyspaurenia too, as well as to prevent stenosis that might occur after sex-transition surgery.



REFs : • A.1004



**Sapigen K** intimate cleanser from Sapi Med has been formulated to deal with problems in gynaecology and proctology. Its surfactants are highly compatible with the skin, and allow delicate yet powerful cleansing of the mucous membranes with good foaming power. The product has strong antibacterial and antimycotics action, and is eudermic and tolerable. The pH is 5.5, which is compatible with the physiological parameters of the intimate membranes.

REFs : • A.3100

sapimed

Products and systems  
for coloproctology

**SAPI MED S.p.A.** is specialized in manufacturing disposable devices for coloproctology and pelvic floor, with more than 30 years experience.

Its product range divides into 3 main lines:

**DIAGNOSTIC**

- Anoscopes
- Rectoscopes
- Ligators

**SURGICAL**

- Surgical anoscopes
- Retractors

**HOME CARE**

- Dilators
- Stomal dilators
- Cleansers
- Creams

# DIAGNOSTIC LINE

sapi med

CE



The **Self•led•light®** disposable proctoscopes are supplied with a built-in LED light, enabling examinations and minor therapeutic procedures to be performed in a non-specialist out-patients department, on the ward or at home. Self•led•light® proctoscopes are made of transparent plastic material and supplied in individually sealed bags. LED lighting can be performed by pulling out the transparent strip, guaranteeing the same light intensity for two hours usage.

REFs

- A.4018L      • A.4022L
- A.4019L      • A.4023L
- A.4021L      • A.4024L

All available in the version without built-in LED light.



Early diagnosis of colorectal cancers and other bowel diseases are made easier with this innovative range of plastic, single-use **rectosigmoscopes** from Sapi Med. Tubes are smooth and made in plastic material with graduation. Each tube has a grip with magnifying lens, connector for insufflator, air insufflator and built-in LED light. By removing the transparent strip the LED powers on, guaranteeing the same light intensity for two hours usage.

REFs

- A.4513LINS    • A.4525LINS
- A.4520LINS

All available in the version without built-in LED light. All available in the STERILE version too.



The LEM **Disposable Haemorrhoid Banding Ligators** are a unique, patented design made from a light, tough polymer. It is a single-patient product available in 4 variants: for use with **suction** unit or **grasping forceps**, **self-suction** and with **4 pre-loaded rings**. The LEMs are designed for use with the Sapi Med Self Light operative proctoscopes.

REFs

- A.5651      • A.5670
- A.5661      • A.5680

# SURGICAL LINE

sapi med

CE



The proctoscopes "**The Beak**" are instruments with a uniquely innovative design feature. The closed, rounded tip allows for introduction of the scope directly into the anus, needing no separate introducer. The proctoscopes "The Beak" offer the specialist the opportunity of clearly and satisfactorily exposing a portion (one fourth of the circumference) of the anal canal and rectal ampulla. They are supplied with built-in LED light. Lighting can be performed by pulling out the transparent strip, guaranteeing the same light intensity for two hours usage.

REFs


- A.4083L      • A.4084L      • A.4084LNS
- All available in the version without built-in LED light.

**ERODE** *EndoRectal Operative Device* is a versatile medical device designed to perform various types of colorectal surgery. In case of procto-pexy, the steps on the plate regulate the degree of needle penetration into the tissue, allowing considerable customisation of the surgical procedure to suit each individual clinical case. The kit includes: an operating anoscope, silk-screened and designed with a window. The opening is 6 cm long with a closed rounded tip and an oval plate with four steps near the tip; one slider with a tab that must be inserted into the anoscope to close completely the window; a knot pusher that can be used when required. It is supplied with built-in LED light. LED Lighting can be performed by pulling out the transparent strip, guaranteeing the same light intensity for two hours usage.



REFs

- A.4087.1L Also available in the version without built-in LED light.

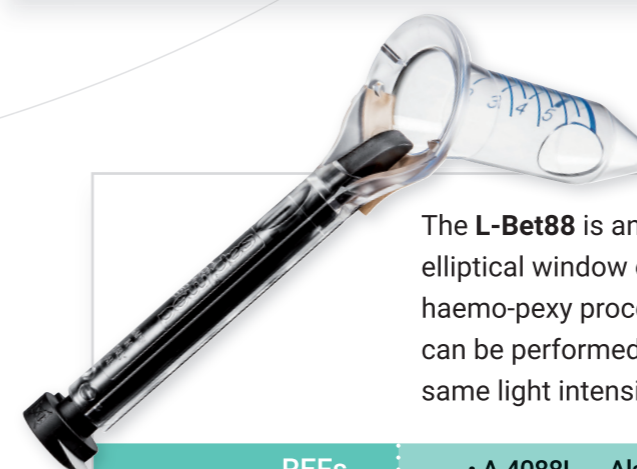

The **L-Bet88** is an operative anoscope featuring a closed tip and an elliptical window on its side. It has been studied to permorm suture haemo-pexy procedure. It is supplied with built-in LED light. LED Lighting can be performed by pulling out the transparent strip, guaranteeing the same light intensity for two hours usage.

REFs

- A.4088L Also available in the version without built-in LED light.