

Contour and straight shape

UDI code and serial number

Low profile design



## Tourniquet Wipe Cuff

Tourniquet cuff, reusable through manual wipe disinfection

# / Tourniquet Wipe Cuff

Tourniquet Wipe Cuff is a reusable tourniquet cuff made of durable material, suitable for repeated use after wipe disinfection. The cuff is available in straight shape (for small or midsize upper extremity) and contour shape (mainly for lower conical shaped extremity) and is provided in different sizes, which are colour coded for optimal cuff selection.

## Features




### Cuff Design

- Low profile design and lightweight
- Large hook and loop fastener for reliable cuff fixation
- Colour coding makes size selection and ordering easier
- Tubing made of kink-resistant silicone material
- Contour and straight shape
- PVC-free



### Special Features

- Specific inner material dries quickly after wipe disinfection and allows timely reuse of the product
- UDI code for product traceability
- Serial number for tracking of product usage
- Processable up to 100 times (wipe disinfection)

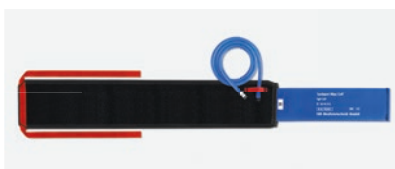


### Straight Shape

- Cylindrical design
- Does not show any difference in diameter between proximal and distal part



### Contour Shape

- Conical design with a smaller diameter distally than proximally when applied
- Occludes the blood flow at lower cuff pressures due to an improved fit to the extremity
- Integrated stiffener to prevent cuff migration

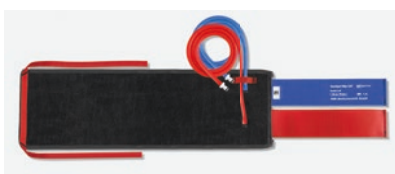

### Double Cuff for Intravenous Regional Anaesthesia (IVRA)

- Double cuff with 2 individually inflatable bladders

## Order information

### Tourniquet Wipe Cuff / Single cuff, reusable

Length	Width	Colour	REF	Box
20 cm	7 cm	□	20-75-700	1
30 cm	8 cm	■	20-75-710	1
35 cm	8 cm	■	20-75-711	1
46 cm	10 cm	■	20-75-712	1
46 cm	12 cm	■	20-75-715*	1
61 cm	12 cm	■	20-75-722*	1
76 cm	12 cm	■	20-75-727*	1
86 cm	12 cm	■	20-75-728*	1
107 cm	12 cm	■	20-75-729*	1

### Tourniquet Wipe Cuff / Double cuff, reusable

Length	Width	Colour	REF	Box
30 cm	12 cm	■	20-77-710	1
46 cm	15 cm	■	20-77-712	1
61 cm	15 cm	■	20-77-722	1

\*Contour shape

This product is manufactured without the use of natural rubber latex, unless otherwise specified.  
This product does not contain phthalates which require labelling according to CLP Regulation (EC) 1272/2008.